

# KOLAY KULLANIM

- Neredeyse acısız<sup>1</sup> delme işlemi sağlayan Clixmotion teknolojisi
- 3 kenarlı delici uçlu özel Softclix lansetleri
- Cilt tipine göre delme derinliği ayarlama seçeneği



## Clixmotion® teknolojisi



- **Doku içinde ve dışında titreşimsiz bir şekilde hareket eder.**
- **Lanset hafifçe durduktan sonra hemen düz bir şekilde geri çekilir.**

### Lanset


3 kenarlı delici uç, farklı cilt tipleri için 11 lanset derinliği ayarı

### Lanset tutucu

Kavrayıcı-tutucu geciktirme mekanizması


### Yönlendirme kanalı

Lansetin titreşimsiz hareketini sağlar

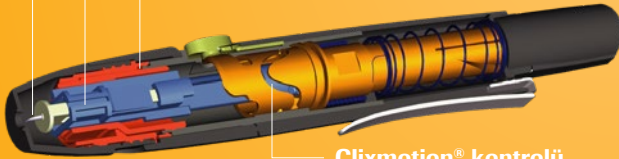

**Clixmotion® kontrolü**  
Soft Stop özelliği ile güvenilir,  
hatasız lanset penetrasyonu

## Diğer lanset teknolojileri

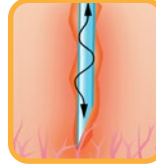

- Lanset cilde girerken titreşir ve ağırlı doku hasarına neden olabilir.
- Lanset aniden durur ve geriye ittirilir.



Başka lansetler, parmak delme cihazının işlevini olumsuz etkileyebilir veya cihaza zarar verebilir ve ağırlı doku hasarına neden olabilir.



Lansetler tek kullanımlıktır.



**Danışma Hattı**  
**0800 211 36 36**

**0850 211 36 36**



[www.diyabetevi.com](http://www.diyabetevi.com)



[facebook.com/diyabetevi](https://facebook.com/diyabetevi)



[instagram/diyabetevi\\_official](https://instagram.com/diyabetevi_official)



[youtube.com/AccuChekTurkey](https://youtube.com/AccuChekTurkey)

### Referans:

1. Kocher S, et al. Comparison of lancing devices for self-monitoring of blood glucose regarding lancing pain. Journal of Diabetes Science and Technology 2009;Volume 3,Issue 5

ACCU-CHEK®, CLIXMOTION® ve SOFTCLIX Roche'a ait ticari markalardır.

© 2021 Roche Diabetes Care